

Evaluation of the Post-Antifungal Effect of Rezafungin and Micafungin against *Candida albicans*, *Candida parapsilosis*, and *Candida glabrata*

Cecilia G. Carvalhaes¹, Jennifer D. Thompson¹, Paul R. Rhomberg¹, Jeffrey B. Locke², Mariana Castanheira¹

¹JMI Laboratories, North Liberty, IA, USA, ²Cidara Therapeutics, San Diego, CA, USA

Introduction

- Post-antifungal effect (PAFE) is defined as the growth suppression of fungal cells after their exposure to an antifungal agent.
- Understanding PAFE helps evaluate dosage regimens for an antifungal agent.
- Rezafungin is a new echinocandin with an extended half-life that exhibits activity against *Candida* spp., *Aspergillus* spp., and *Pneumocystis* spp. and is in Phase 3 clinical development.
- In this study, the PAFE of rezafungin was compared to that of micafungin against *Candida albicans*, *C. glabrata*, and *C. parapsilosis* isolates.

Materials and Methods

- Six *Candida* spp. isolates were tested, including 2 *C. albicans* (ATCC 90028 and #1 clinical isolate), 2 *C. parapsilosis* (ATCC 22019 and #2 clinical isolate), and 2 *C. glabrata* (#3 and #4 clinical isolates).
- Antifungal susceptibility testing was performed using the Clinical and Laboratory Standards Institute (CLSI) reference broth microdilution method.
 - Rezafungin (Cidara Therapeutics) and micafungin (Sigma-Aldrich) were tested.
 - Each isolate was tested in triplicate to establish baseline MIC values (bMIC).
 - Modal MIC values were used to determine drug concentrations for PAFE testing.
- For PAFE determinations, antifungal concentrations of 1X, 4X, and 16X the baseline MIC were used.
- A starting inoculum of 1-5x10⁵ CFU/mL from a fresh culture was added to RPMI with the respective antifungal at the desired concentration (1X, 4X, or 16X the bMIC).
 - A growth control that was not exposed to an antifungal agent was used to determine the standard 1-log₁₀ increase.
- After a 1-h exposure to the antifungal agent, the cells were washed three times with RPMI and reconstituted with pre-warmed RPMI to a final volume of 10mL.
- Colony counts were performed at T0 (pre-exposure), after the 1-h drug exposure, and after the cell wash (T1).
 - Test cultures were re-incubated following final cell wash.
 - Colony counts were performed at T2, T4, T8, T12, T24, and T48 hours.
- PAFE was calculated as the difference in time required for isolates to grow 1-log₁₀ after the final cell wash compared to the untreated growth control.
 - PAFE = T – C, where T is time required for isolate to increase 1-log₁₀ (in CFU) after drug removal and C is time required for untreated growth control to increase 1-log₁₀ after undergoing the same process performed on the test culture.
- The reduction in starting inocula, in log-kill, was also calculated over the 48-h study period.

Results

- PAFE and bMIC results for rezafungin and micafungin are shown in Table 1.
- *C. albicans*
 - Rezafungin and micafungin PAFEs were >14.9 h against the *C. albicans* clinical isolate for all concentrations tested.
 - The *C. albicans* ATCC 90028 control failed to re-grow over 1-log₁₀; the rezafungin and micafungin PAFEs could not be determined against this isolate.
- *C. glabrata*
 - Rezafungin PAFE results were >40 h for both *C. glabrata* strains, regardless of the concentration tested.
 - Micafungin PAFE results were equivalent to rezafungin PAFE values (>40 h) for *C. glabrata* isolates at all concentrations, except the 1X bMIC for the *C. glabrata* #3 isolate.
 - The micafungin PAFE value was 20.4 h for *C. glabrata* #3 isolate, while the rezafungin PAFE was >46.7 h.

C. parapsilosis

- Rezafungin PAFE results were also >40 h against the *C. parapsilosis* ATCC 22019 strain at all concentrations.
- In contrast, no micafungin PAFE was observed against *C. parapsilosis* ATCC 22019 at 1X and 4X bMIC, and a short PAFE (1.8 h) was noted at 16X bMIC.
- The *C. parapsilosis* clinical isolate #2 displayed prolonged rezafungin PAFE values (range, 18.4 h to >36.6 h), regardless of the concentration tested.
- A short PAFE was displayed by micafungin at 1X (1.6 h) and 4X (7.4 h) bMIC against *C. parapsilosis* #3 while a PAFE of 31.3 h was noted for micafungin against this clinical isolate at 16X bMIC.

Conclusion

- Rezafungin showed sustained growth inhibition following drug removal and displayed equivalent or longer PAFE values than micafungin against all tested *Candida* spp.

Acknowledgements

This study was supported by Cidara Therapeutics.

References

CLSI. M27Ed4E. Reference Method for Broth Microdilution Antifungal Susceptibility Testing of Yeasts, fourth edition. Wayne, PA, Clinical and Laboratory Standards Institute, 2017.

CLSI. M60Ed2E. Performance standards for antifungal susceptibility testing of yeasts. Wayne, PA, Clinical and Laboratory Standards Institute, 2020.

Ernst EJ, Klepser ME, Pfaller MA. Postantifungal effects of echinocandin, azole, and polyene antifungal agents against *Candida albicans* and *Cryptococcus neoformans*. *Antimicrob Agents Chemother.* 2000;44(4):1108–1111.

Nguyen KT, Ta P, Hoang BT, et al. Characterising the post-antifungal effects of micafungin against *Candida albicans*, *Candida glabrata*, *Candida parapsilosis* and *Candida krusei* isolates. *Int J Antimicrob Agents.* 2010;35(1):80–84.


Shields RK, Nguyen MH, Press EG, Clancy CJ. Five-minute exposure to caspofungin results in prolonged postantifungal effects and eliminates the paradoxical growth of *Candida albicans*. *Antimicrob Agents Chemother.* 2011 Jul;55(7):3598-602.

Smith RP, Baltch A, Bopp LH, Ritz WJ, Michelsen PR. Post-antifungal effects and time-kill studies of anidulafungin, caspofungin, and micafungin against *Candida glabrata* and *Candida parapsilosis*. *Diagn Microbiol Infect Dis.* 2011;71(2):131–138.

Contact


Cecilia Carvalhaes, MD, PhD
 JMI Laboratories
 345 Beaver Kreek Centre, Suite A
 North Liberty, Iowa 52317
 Phone: (319) 665-3370
 Fax: (319) 665-3371
 Email: Cecilia-carvalhaes@jmilabs.com

Figure 1. Rezafungin PAFE against *C. albicans*, *C. parapsilosis*, and *C. glabrata* at 1X, 4X, and 16X the baseline MIC values

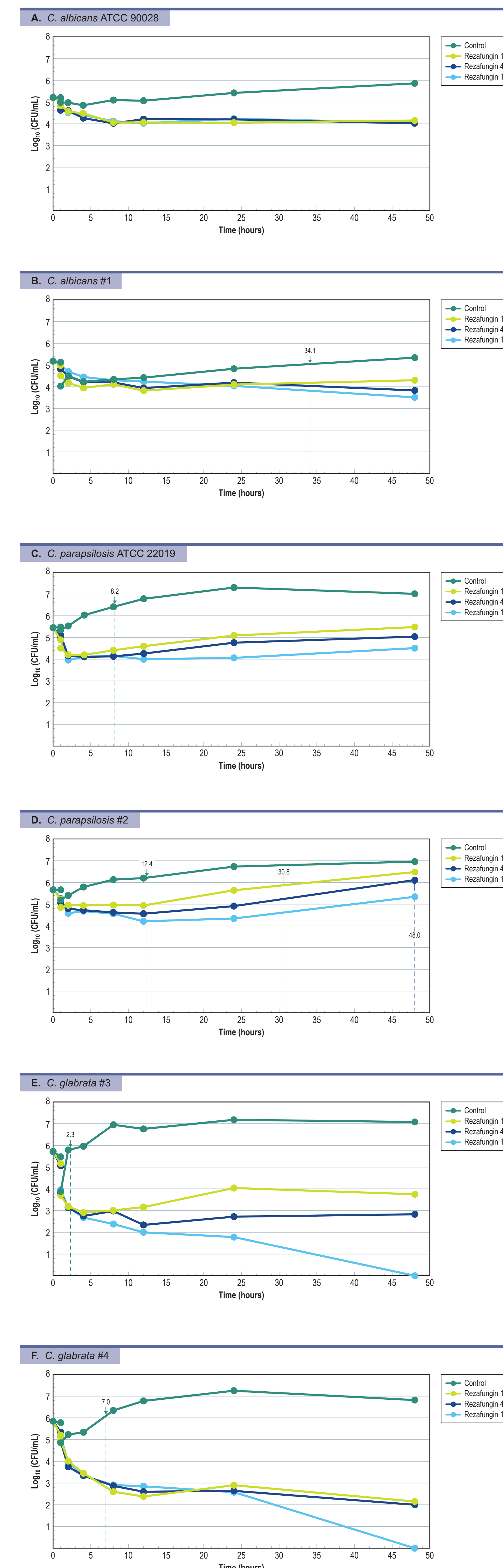

Dashed lines indicate the time required for the isolate to increase 1-log₁₀ (in CFU) after drug removal

Figure 2. Micafungin PAFE against *C. albicans*, *C. glabrata*, and *C. parapsilosis* at 1X, 4X, and 16X the baseline MIC values

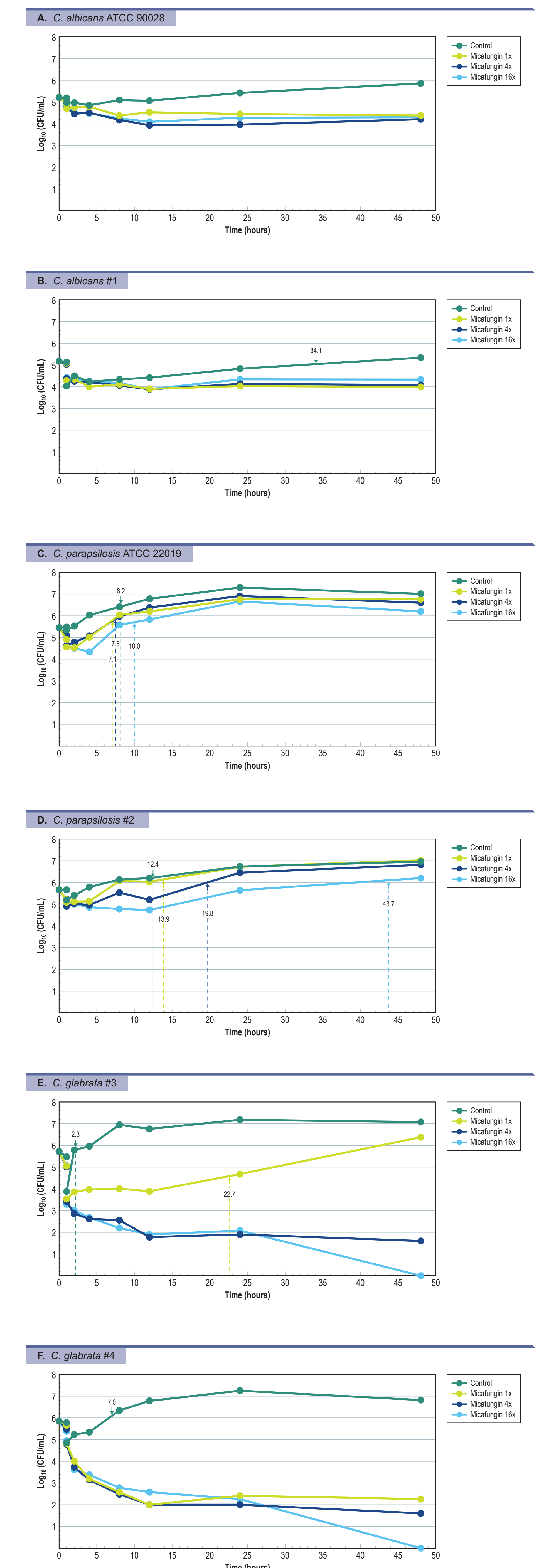

Dashed lines indicate the time required for the isolate to increase 1-log₁₀ (in CFU) after drug removal

Table 1. PAFE for rezafungin and micafungin against *C. albicans*, *C. parapsilosis*, and *C. glabrata*

Antifungal/Strain	Baseline MIC (mg/L)	PAFE (hours) at the following multiple of baseline MIC		
		1X	4X	16X
Rezafungin				
<i>C. albicans</i> ATCC 90028	0.03	ND	ND	ND
<i>C. albicans</i> #1	0.06	>14.9	>14.9	>14.9
<i>C. parapsilosis</i> ATCC 22019	1	>40.8	>40.8	>40.8
<i>C. parapsilosis</i> #2	1	18.4	35.6	>36.6
<i>C. glabrata</i> #3	0.12	>46.7	>46.7	>46.7
<i>C. glabrata</i> #4	0.12	>42.0	>42.0	>42.0
Micafungin				
<i>C. albicans</i> ATCC 90028	0.03	ND	ND	ND
<i>C. albicans</i> #1	0.015	>14.9	>14.9	>14.9
<i>C. parapsilosis</i> ATCC 22019	1	≤0.0	≤0.0	1.8
<i>C. parapsilosis</i> #2	1	1.6	7.4	31.3
<i>C. glabrata</i> #3	0.06	20.4	>46.7	>46.7
<i>C. glabrata</i> #4	0.03	>42.0	>42.0	>42.0

Abbreviations: ND, not determined